KLOPFY

NM-088

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ENGLISH

The ORIGINAL

User manual for
tapping & shiatsu massager
KLOPFY NM-088
WITH ZIPPER

Z2 091698 0009 Rev.03

KLOPFY

Contents:

Introduction	15
Explanation - signs and symbols	15
Traditional shiatsu and tapping massage	16
Signs and symbols	16
Intended use	16
Getting to know your device	16
Parts of the product	17
Product functions	18
Massage positions	18,10
Safety precautions	19
Specific group of people	20
Contents included	21
Technical specifications	21
Maintenance and storage	22
Electric and electronic devices	22
Regulatory compliance	23
Warranty	24

Introduction

Dear customer,

Thank you for choosing Donnerber Original KLOPFY neck and shoulder massager. Donnerberg Original KLOPFY massager is manufactured to the highest German quality standards certified by TÜV.

Please read the manual instructions carefully and save it for later use, be sure to make them accessible to the other users. Please follow the instructions and use the product as described.

Donnerberg Team

Explanation - SIGNS AND SYMBOLS

WARNING – Warning instruction indicating a risk of injury

IMPORTANT – Safety note regarding potential for damage to the device

Traditional shiatsu and tapping massage

Tapping massage & shiatsu massage combined in one device. Tapping is a light percussion (also known as tapotement) movement.

The pleasant and penetrating tapping massage relaxes the muscles and relieves tension. The word Shiatsu means 'finger pressure' in Japanese. Shiatsu uses pressure, touch and manipulative techniques to bring the body and its inner energy paths into balance. Shiatsu massage can be used to reduce stress and promote overall health and well-being.

Intended use

The unit is only intended for domestic /private use, not for commercial use. It is designed for massaging the following body parts: neck, back, shoulders, belly, legs and arms. If you are unsure whether the neck and shoulder massager is suitable for you, please consult your doctor.


Getting to know your device

USAGE AND DIRECTION - Sit in a comfortable chair, carefully place the massager over your shoulders and switch it on. Put the arms through the loops of arm straps to facilitate the positioning of the massager and to control the intensity of the massage. Place the massager on the desired areas or parts of the body and adjust the intensity of the massage by pulling harder or softer on the loops.

Each massage cycle lasts 15 minutes but can be shortened on demand. When the 15 minute cycle has elapsed, the massager will automatically switch off. This device must be turned off when not in use.


Parts of the product

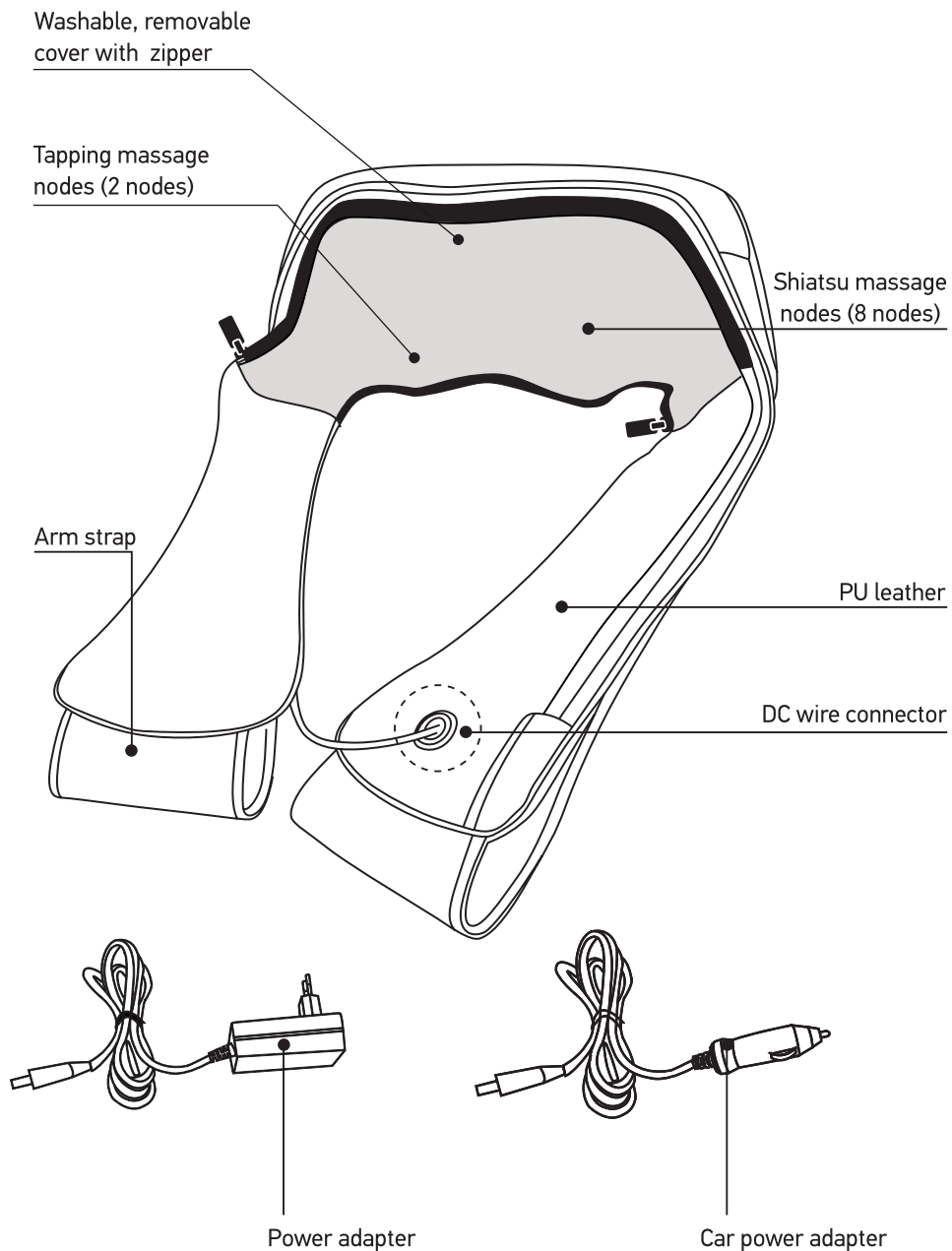Product functions

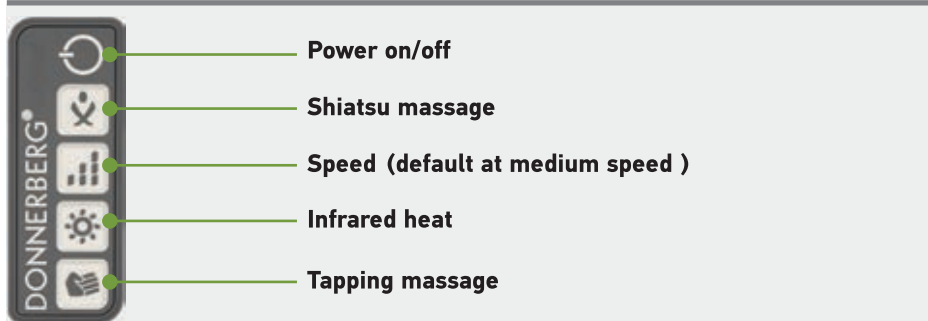Power on/off

- Press this button once to turn on the device
- The medium massage speed and high intensity infrared heat is on by default.
- One minute shiatsu, one minute shiatsu combined with tapping massage will perform and each function will shift after one minute.

Shiatsu massage

Press this button to lock - start shiatsu massage. This function is valid for shiatsu massage function only, not for tapping massage.

Speed (default at medium speed)

- Press this button once for low speed;
- Press this button twice for high speed;
- Cycles from medium to low to high speed will continuously repeat in this sequence.

Infrared heat

Infrared heat only

Heating function can be controlled independently without massage

- High-intensity Infrared heat turns on when you press this button once
- Press this button twice to turn on low-intensity Infrared heat
- Press this button three times to turn off Infrared heat

Combination of Massage & Infrared heat

- Press the Power On/Off button to turn on default mode (interchangeably perform shiatsu function or shiatsu&tapping function each minute with high intensity Infrared heat)
- Press Infrared heat button once to change heat intensity from high to low
- Press Infrared heat button twice to turn off this function

Tapping massage


Press the button to lock - start shiatsu & tapping function. The massager will perform tapping & shiatsu all the time.
The massager will be switched off automatically after 15 minutes

Massage positions

Please apply on body parts shown in the pictures below

Neck & Shoulders

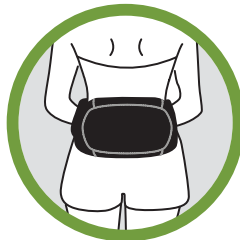
Back

Belly

Legs

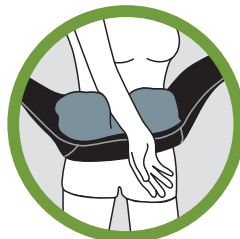Arms

Safety precautions

To prevent the risk of fire, electric shock, or other injury, please follow the next points:

1. Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plates is compatible with your mains supply. Only insert the power supply unit into a socket when the appliance is switched off.
2. Do not cover or use the unit under a blanket when operating. Do not use while bathing or showering. The unit must be kept dry and stored in a dry and safe place. The power socket,

cable, and appliance must be at appropriate distance from the heat, hot surfaces, moistures and liquids. Connect the unit only in that way, so that the power plug is freely accessible. The power plug should be freely accessible when inserted into the socket. The cable must not be kinked, trapped or twisted.

3. To use this product safely and correctly, read this manual thoroughly and operate the device as instructed. Never operate this appliance if it has a damaged cord or plug or if it is not working properly. In case the device is not working or has a fault, send it to the service centre for examination. Do not plug in the device in the damaged sockets. If the device has been dropped into water, do not pick it up or touch it. Pull out the power plug from the socket immediately.

Do not operate where aerosol (spray) products are being used or where oxygen is being administered. After you have finished using the unit turn off the power. Before unplugging, the device must be turned off. Do not obstruct with objects or thrust.

Do not use the massager while driving!

Children are not allowed to play with the device.

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation and choking hazard!

SPECIFIC GROUP OF PEOPLE

- This appliance can be used by children aged 8 and over and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved.
- Do not use this device to supplement or replace medical treatment
- Chronic complaints and symptoms may be made worse.

You should refrain from using the massager, or should consult your doctor before using it, if

- You are pregnant,
 - You have a cardiac pacemaker, artificial joints or electronic implants,
 - You suffer from one or more of the following illnesses and complaints: Circulation disorders, varicose veins, open wounds, bruises, broken skin, inflammation of the veins.
-
- Do not use the appliance near your eyes or near any other sensitive areas.
 - Take care when using the heating function. Do not use if you have circulatory problems or on areas sensitive to heat.
 - The surface of the appliance is warm. People who are sensitive to heat must be careful when using the appliance.
 - If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
 - Before using the massager, consult your doctor regarding any unexplained pain, if you are undergoing medical treatment or using medical equipment

Contents included

1x Donnerberg Original Klopfy neck and shoulder massager NM-088
1x Power adapter (length approx 2.8m)
1x Car charger
2x Removable cover
1x Soft cloth cover for sensitive skin
1x Protective custody - environmentally friendly
1x User manual

Technical specifications


Product name	KLOPFY NM - 088
Product model	NM-088 Black
Voltage supply	Output DC12V 3A 36W
Input voltage	100-240 V, 50/60 Hz
Automatic switch off	After approx. 15 minutes
Net Weight	Approx. 2,2 kg
Gross Weight	Approx. 2,96 kg
Dimensions	53 x 42 x 20,5 cm

Maintenance and storage

- Before cleaning, make sure the switch is turned off and power plug is disconnected from the outlet. Do not use abrasive cleaners such as benzene or thinner to clean the product, otherwise the colour of the product may fade. Do not immerse the appliance in the water or the other liquids during the cleaning! Small marks can be removed by using a damp cloth. Please remember to clean only the top side of the device.
- You may only clean the massager yourself. If the device malfunctions please do not repair it yourself as this leads to automatic expiry of all guarantee rights.
- Thanks to the integrated zipper, the cover can be easily removed and washed in the washing machine at 30° C. Simple and easy, just unzip the zipper and replace the cover with another one.

We would kindly advise you to clean your massager regularly by taking off the cover and removing dirt (hair, dust) that gathered between the massage node and the contact plate with cotton buds. Regular removal of contamination can prevent malfunctions of certain components of infrared heat and massage.

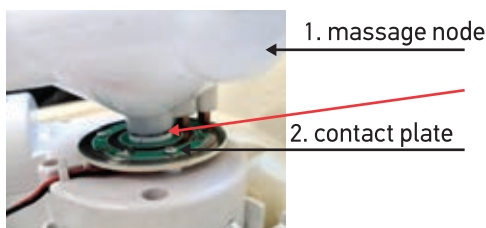Electric and electronic devices

Information for private households

The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE) in the United Kingdom contains several requirements concerning the handling of electric and electronic devices as waste. In the following paragraphs, we list the most important requirements for you.

1. Separated recycling of „old devices“

Electric and electronic devices that have become waste are called old devices. Owners of such old devices have to bring these to recycling facilities, where these are handled and recorded separately from urban solid waste. Old devices must not end up in the household refuse. In fact they have to be collected in special collecting and recycling systems.

2. Possibility for disposal of „old devices“

Owners of old devices from private households can dispose these at the public waste disposal authorities or bring them to producer or seller related disposal locations. You can find these disposal locations online: <https://www.recycle-more.co.uk/>

3. Privacy information

Old devices often contain private data, especially appliances from data processing and communication equipment, as smartphones or hard discs. In your own interest, please arrange to eliminate all personal and private data from your disposed equipment.

4. Meaning of the symbol „crossed out waste container“

The labelled symbol showing the „crossed out waste container“ on electric and electronic equipment indicates that

the concerning equipment has to be disposed separately from urban solid waste.

5. The Waste Electrical and Electronic Equipment Regulations 2018 (WEEE)

Following the implementation of the Directive 2012/19/EU and its British transportation WEEE Regulations 2018, we as a producer of EEE are registered in the Public Register of Electrical and Electronic Equipment Producers. The register is held in England and Wales by the Environment Agency.

Regulatory compliance

CE (Conformité Européenne)- indicates that the product conforms with relevant EU directives regarding health and safety or environmental protection

TÜV (Technical Inspection Association)-product was tested for electrical safety, hazardous material components, and unwanted noise production.

Z2 091698 0009 Rev. 02

RoHS (Restriction of Hazardous Substances)- restricts the use of specific hazardous materials found in electrical and electronic products.

Warranty

The company Donnerberg, Mr. Svetozar Grbic, Von-Werthern-Str. 4a, 96487 Dörfles-Esbach, Germany grants a valid manufacturer's warranty for the territory of the European Union and the civil-law defect liability. The warranty starts from the time the device was delivered to the customer. It is not extended by a previous warranty case. The warranty is granted only to the buyer and is not transferable to third parties.

In addition to the statutory liability for defects, we also repair or replace equipment whose defects are not initial defects unless they are caused by improper use or are insignificant. The choice between replacement and repair is ours. The warranty is only granted if the defect is notified immediately upon its occurrence and is reported. You must include the claim form and the copy of your purchase receipt as the proof of warranty. Please note that we will not process warranty claims without this claim form.

You may report your warranty claim to one of the following:

Company Donnerberg
Von-Werthern-Str. 4a, 96487 Dörfles-Esbach
Germany

de.info@donnerberg.net
or telephone: + 49 89 215 400 980 / + 49 95 619 769 199

Upon your notification, our customer service team will send you a claim form.

For questions regarding warranty or any other questions, you may also contact us via email:
uk.info@donnerberg.net

Donnerberg S.G.

Von-Werthern-Str. 4a

96487 Dörfles-Esbach, Deutschland

Tel: + 49 89 215 400 980 / + 49 95 619 769 199

Mail: de.info@donnerberg.net

Web: www.donnerberg.net